

NARCOBIT-E (Type II)

Popular



KN-1071-I
For isoflurane



KN-1071-S
For sevoflurane

- **Animals for use**

The device is designed for use mainly in small experimental animals such as mice and rats.

- **Anesthetics to be used**

Anesthetic can be selected from isoflurane and sevoflurane. *The anesthetic device is designed for use with one of these drugs exclusively.

- **Simultaneous maintenance of multiple animals**

By combining peripheral apparatuses, 1 to 10 mice or 1 to 8 rats can be maintained.

- **Necessary functions, only as required**

We have achieved further cost reductions compared to our previous products by reviewing each component. We provide a wide variety of peripheral apparatuses to answer various needs.

Main specifications	Main body size : W 175 x D 145 x H 360 mm
	Main body weight : Approx. 7.7 kg (excluding the Colorful-o-meter)
	Vaporizer : For isoflurane or sevoflurane (for use with one anesthetic drug exclusively)
	Flow rate : Induction anesthesia box side: Approx. 5 L/min*1 Maintenance anesthesia mask side: Up to 1.0 L/min*2
	Power source : AC 100 V

Accessories³ KN-1013-1 Colorful-o-meter 1 L type x1
(For the control of flow rate in maintenance anesthesia mask)



*1: Induction anesthesia box can be filled with anesthetic gas.


If you wish to adjust the flow rate, we have KN-1013-5 Colorful-o-meter (5 L type) x 1 (sold separately) available to be connected to the induction anesthesia box side.

*2: One Colorful-o-meter is provided with the device for use on the side of maintenance anesthetic mask.

For the maintenance anesthesia in multiple animals, Colorful-o-meter needs to be provided for each animal to be maintained.

*3: In addition to this, induction anesthesia box and other apparatus will be required (sold separately). We have a catalogue available for "related products for inhalation anesthesia apparatus".

Connection of related products for KN-1071 NARCOBIT-E (Type II)

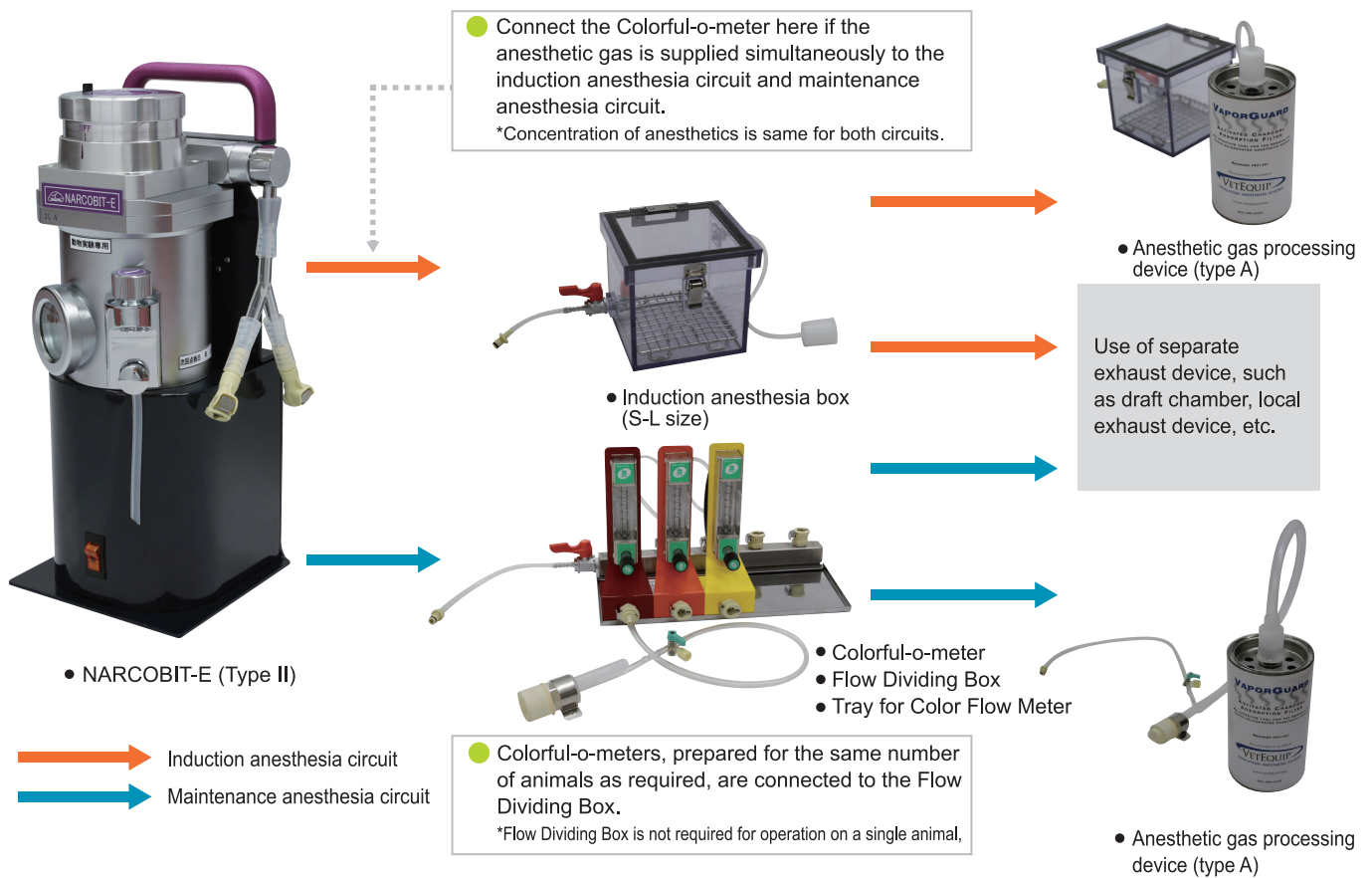

■ Example of set (for single animal)

Products used

1. KN-1071 NARCOBIT-E (Type II)
2. KN-1010 Induction anesthesia box set
3. KN-1013 Colorful-o-meter
4. KN-1019 Natsume type anesthetic mask
5. KN-1050 Anesthetic gas processing device (type A)

*Please inquire for more details, as we have a catalog available for related products for inhalation anesthesia apparatus, as well as catalogues for individual products.

■ Example of set (for multiple animals)

Products used

1. KN-1071 NARCOBIT-E (Type II)
2. KN-1010 Induction anesthesia box set
3. KN-1013 Colorful-o-meter
4. KN-1013-B5 Flow Dividing Box
5. KN-1013-T5 Tray for Color Flow Meter
6. KN-1013-ST Stand for containing the anesthesia box
7. KN-1019 Natsume type anesthetic mask
8. KN-1050 Anesthetic gas processing device (type A)

*Please be aware that appearance and specifications may change for improvement without a preliminary announcement.

Please feel free to inquire us:
Please let us know the animals to be used (mice? rats?), number of animals to be anesthetized simultaneously, and whether induction anesthesia and maintenance anesthesia will be used simultaneously.
We will suggest the anesthetics device for you, including the related products.



株式会社 夏目製作所

Headquarters 2-18-6 Yushima Bunkyo-ku, Tokyo 113-8551
Tel: 03-3813-3251 Fax: 03-3815-2002
HP: <https://www.nazme.co.jp>

Physics and chemistry apparatus / Basic medical apparatus / Pharmaceutical research apparatus
Laboratory animal breeding management apparatus / Medical science apparatus general

